



**ohmic instruments co.**

508 August Street, Easton, Maryland 21601  
(410) 820-5111 (800) 626-7713 Fax (410) 822-9633  
www.ohmicinstruments.com  
e-mail: ohmic@ohmicinstruments.com



**BIOMEDICAL & ENVIRONMENTAL SENSORS, INSTRUMENTS & CONTROLS SINCE 1969**


**MULTIPURPOSE HOLDER - POSITIONER MODEL MHP-102  
REPLACES PREVIOUS TRANSDUCER HOLDER FOR ALL UPM-DT SERIES**

**THE NEW! Model: MHP-102  
Multipurpose Holder Positioner**

- SCREWS TO BASE IN PLACE OF CLAMP ASSEMBLY
- GRIP RANGE: TO 2 1/2 INCHES
- LOCKING BALL JOINTS EASILY ADJUST
- LARGE THUMB SCREWS ON CLAMP



PREVIOUS CLAMP ASSEMBLY

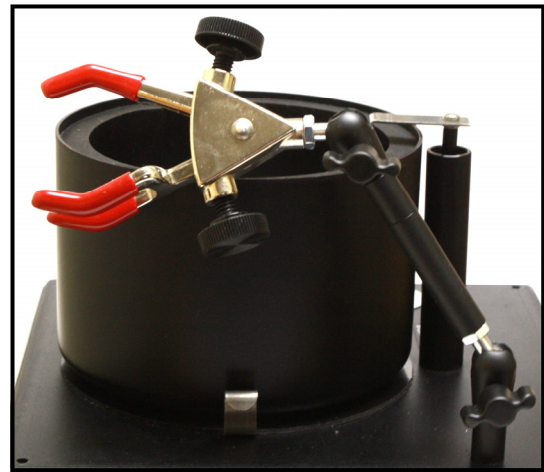

NEW CLAMP! MHP-102 INSTALLED

**EASY FAST ACCURATE POSITIONING OF TRANSDUCERS, PROBES AND TEST EQUIPMENT**

OHMIC Instruments' Model MHP-102 is the current clamping/holding device used on all OHMIC UPM-DT's. The MHP-102 consists of a sturdy clamp head coupled to a locking ball joint, extension shaft, and a second locking ball joint. THE MHP-102 mounts directly to the platform and replaces the clamp assembly and support rod on all OHMIC UPM-DT series Ultrasound Power Meters.



Also available is the MHP-101 block clamp assembly. The anodized block adapts the new clamp to the existing target support rod.



"MHP-101 Installed"

The ball joint positioner can be locked down or adjusted to allow pinpoint adjustments to be made. Both the MHP-101 and MHP-102 provide easy replacement and upgrading of the clamp assembly on your OHMIC Ultrasound Power Meter.

**ENGINEERING SUPPORT:** OHMIC Instruments Co. designs and manufactures a full line of environmental and bio-medical sensors, instruments and controls. Many of our products are custom designed to meet specific requirements. Our engineers will be pleased to discuss your application.

**WARRANTY:** All products manufactured by OHMIC Co. are warranted to be free of defects in material and workmanship for one year after delivery. Any equipment found to be defective within this period will be repaired or replaced free of charge.